



A Sismex Group Company



### Instructions For Use

REF: LPH012-S / LPH012 / LPH012-20

## TEL/AML1 (ETV6/RUNX1) Translocation, Dual Fusion Probe



PROFESSIONAL USE ONLY



www.cytocell.com

Further information and other languages available at [www.cytocell.com](http://www.cytocell.com)

### Limitations

This device is designed to detect rearrangements with breakpoints in the region covered by the red and green clones in this probe set, which includes the *TEL* (*ETV6*) and *AML1* (*RUNX1*) regions. Breakpoints outside this region, or variant rearrangements wholly contained within this region, may not be detected with this product.

The test is not intended for: use as a stand-alone diagnostic, prenatal testing, population-based screening, near-patient testing or self-testing. This product is intended for laboratory professional use only; all results should be interpreted by suitably-qualified staff, taking into account other relevant test results.

This product has not been validated for use on sample types or disease types other than those specified in the intended use.

Reporting and interpretation of FISH results should be consistent with professional standards of practice and should take into consideration other clinical and diagnostic information. This kit is intended as an adjunct to other diagnostic laboratory tests and therapeutic action should not be initiated on the basis of the FISH result alone.

Failure to adhere to the protocol may affect the performance and lead to false positive/negative results.

This kit has not been validated for purposes outside of the intended use stated.

### Intended Use

The Cytocell® Aquarius TEL/AML1 (*ETV6/RUNX1*) Translocation, Dual Fusion Probe is a qualitative, non-automated, fluorescence *in situ* hybridisation (FISH) test used to detect chromosomal rearrangements between the 12p13.2 region on chromosome 12 and the 21q22.1 region on chromosome 21 in Carnoy's solution (3:1 methanol/acetic acid) fixed haematologically-derived cell suspensions from patients with confirmed or suspected acute lymphoblastic leukaemia (ALL).

### Indications

This product is designed as an adjunct to other clinical and histopathological tests in recognised diagnostic and clinical care pathways, where knowledge of *TEL-AML1* (*ETV6-RUNX1*) translocation status would be important for clinical management.

### Principles of the Test

Fluorescence *in situ* hybridisation (FISH) is a technique that allows DNA sequences to be detected on metaphase chromosomes or in interphase nuclei from fixed cytogenetic samples. The technique uses DNA probes that hybridise to entire chromosomes or single unique sequences, and serves as a powerful adjunct to G-banded cytogenetic analysis. This technique can now be applied as an essential investigative tool within prenatal, haematological and solid tumour chromosomal analysis. Target DNA, after fixation and denaturation, is available for annealing to a similarly denatured, fluorescently labelled DNA probe, which has a complementary sequence. Following hybridisation, unbound and non-specifically bound DNA probe is removed and the DNA is counterstained for visualisation. Fluorescence microscopy then allows the visualisation of the hybridised probe on the target material.

### Probe Information

The cytogenetically-cryptic t(12;21)(p13;q22) translocation between *ETV6* (*ets variant 6*) at 12p13 and *RUNX1*, (*RUNX family transcription factor 1*) at 21q22, results in the *ETV6-RUNX1* chimeric fusion gene<sup>1</sup>.

The *ETV6* and *RUNX1* genes both encode transcription factors; *ETV6* has been shown to be required for proper transcription during haematopoiesis within the bone marrow<sup>1,2</sup>. The *ETV6-RUNX1* protein converts *RUNX1* to a transcriptional repressor and causes overexpression of the erythropoietin receptor (EPOR) and activation of downstream JAK-STAT signaling<sup>1</sup>.


B-lymphoblastic leukaemia/lymphomas with t(12;21)(p13;q22) translocations form a recognised disease entity according to the World Health Organization (WHO) classification of myeloid neoplasms and acute leukaemia. This is the most common sub-group of childhood B-ALL accounting for about 25% of cases<sup>3</sup>. As the t(12;21)(p13;q22) translocation is cytogenetically-cryptic, FISH is an important diagnostic tool for this leukaemia<sup>4</sup>.

B-ALL with *ETV6-RUNX1* is considered to have a favourable outcome with cure rates more than 90%<sup>3</sup>. Late relapses have been reported; these have been attributed to the presence of persistent preleukaemic clones that survived chemotherapy<sup>3,5</sup>.

*ETV6* has also been shown to be deleted in some children with ALL, with loss of heterozygosity (LOH) of chromosome 12p12-13; these deletions often seen in the presence of *ETV6-RUNX1* translocations<sup>6</sup>.

### Probe Specification

TEL, 12p13.2, Red  
AML1, 21q22.12, Green

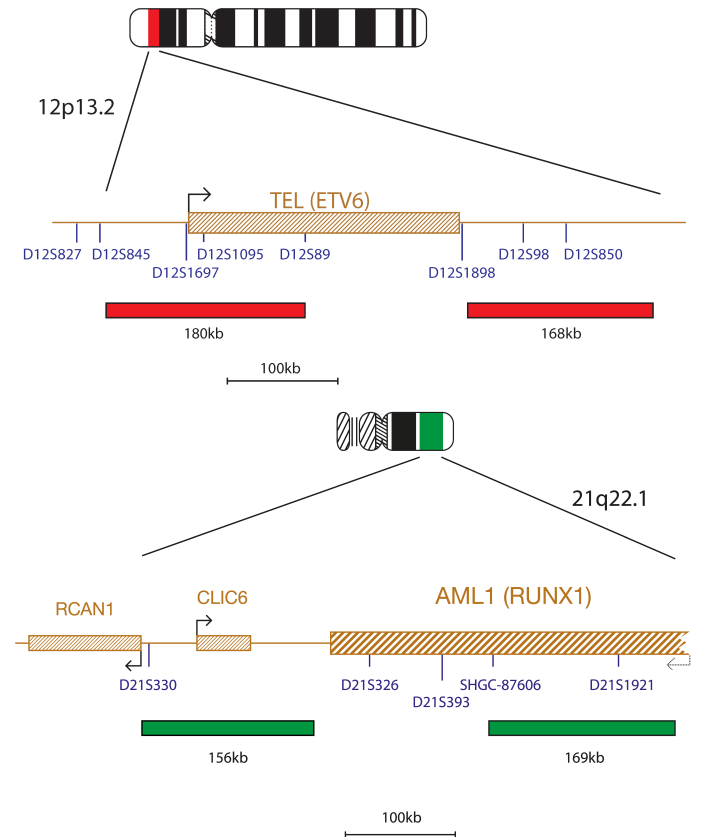

The TEL probe mix, labelled in red, contains a probe covering a 180kb region between the markers D12S845 and D12S89 and a second probe centromeric to the TEL (*ETV6*) gene, extending 168kb from the marker D12S1898. For AML1, there are two probes, labelled in green, one covering a 156kb region centromeric to the AML1 (*RUNX1*) gene, including the *CLIC6* gene and a second probe covering a 169kb region, including the markers D21S1895 and D21S1921.

### Materials Provided

**Probe:** 50µl per vial (5 tests), 100µl per vial (10 tests) or 200µl per vial (20 tests)  
The probes are provided premixed in hybridisation solution (formamide; dextran sulphate; saline-sodium citrate (SSC)) and are ready to use.

**Counterstain:** 150µl per vial (15 tests) or 500µl per vial (50 tests)


The counterstain is DAPI antifade (ES: 0.125µg/ml DAPI (4,6-diamidino-2-phenylindole)).

### Warnings and Precautions

1. For *in vitro* diagnostic use. For professional use only.
2. Wear gloves when handling DNA probes and DAPI counterstain.
3. Probe mixtures contain formamide, which is a teratogen; do not breathe fumes or allow skin contact. Wear gloves, a lab coat, and handle in a fume hood. Upon disposal, flush with a large volume of water.

- DAPI is a potential carcinogen. Handle with care; wear gloves and a lab coat. Upon disposal, flush with a large volume of water.
- Dispose of all hazardous materials according to your institution's guidelines for hazardous waste disposal.
- Operators must be capable of distinguishing the colours red, blue and green.
- Failure to adhere to the outlined protocol and reagents may affect the performance and lead to false positive/negative results.
- The probe should not be diluted or mixed with other probes.
- Failure to use 10µl of probe during the pre-denaturation stage of the protocol may affect the performance and lead to false positive/negative results.

#### Storage and Handling

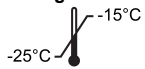

The Aquarius® kit should be stored between -25°C to -15°C in a freezer until the expiry date indicated on the kit label. The probe and counterstain vials must be stored in the dark.



The probe remains stable throughout the freeze-thaw cycles experienced during normal use (where one cycle constitutes the probe's removal from and replacement into the freezer) and is photostable for up to 48 hours after being exposed to continuous lighting conditions. All efforts must be made to limit exposure to light and temperature changes.

#### Equipment and Materials Necessary but not Supplied

Calibrated equipment must be used:

- Hotplate (with a solid plate and accurate temperature control up to 80°C)
- Calibrated variable volume micropipettes and tips range 1µl - 200µl
- Water bath with accurate temperature control at 37°C and 72°C
- Microcentrifuge tubes (0.5ml)
- Fluorescence microscope (Please see Fluorescence Microscope Recommendation section)
- Phase contrast microscope
- Clean plastic, ceramic or heat-resistant glass Coplin jars
- Forceps
- Calibrated pH meter (or pH indicator strips capable of measuring pH 6.5 – 8.0)
- Humidified container
- Fluorescence grade microscope lens immersion oil
- Bench top centrifuge
- Microscope slides
- 24x24mm coverslips
- Timer
- 37°C incubator
- Rubber solution glue
- Vortex mixer
- Graduated cylinders
- Magnetic stirrer
- Calibrated thermometer

#### Optional Equipment not Supplied

- Cytogenetic drying chamber

#### Reagents Needed but not Supplied

- 20x saline-sodium citrate (SSC) Solution
- 100% Ethanol
- Tween-20
- 1M Sodium hydroxide (NaOH)
- 1M Hydrochloric acid (HCl)
- Purified water

#### Fluorescence Microscope Recommendation

Use a 100-watt mercury lamp or equivalent and oil immersion plan apochromat objectives 60/63x or 100x for optimal visualisation. The fluorophores used in this probe set will excite and emit at the following wavelengths:

Fluorophore	Excitation <sub>max</sub> [nm]	Emission <sub>max</sub> [nm]
Green	495	521
Red	596	615

Ensure appropriate excitation and emission filters that cover the wavelengths listed above are fitted to the microscope. Use a triple bandpass DAPI/green spectrum/red spectrum filter or a dual bandpass green spectrum/red spectrum filter for optimal simultaneous visualisation of the green and red fluorophores.

Check the fluorescence microscope before use to ensure it is operating correctly. Use immersion oil that is suitable for fluorescence microscopy and formulated for low auto fluorescence. Avoid mixing DAPI antifade with microscope immersion oil as this will obscure signals. Follow manufacturers' recommendations in regards to the life of the lamp and the age of the filters.

#### Sample Preparation

The kit is designed for use on haematologically-derived cell suspensions fixed in Carnoy's solution (3:1 methanol/acetic acid) fixative, that are prepared according to the laboratory or institution guidelines. Prepare air dried samples on microscope slides according to standard cytogenetic procedures. The AGT *Cytogenetics Laboratory Manual* contains recommendations for specimen collection, culturing, harvesting and for slide making<sup>6</sup>.

#### Solution Preparation

##### Ethanol Solutions

Dilute 100% ethanol with purified water using the following ratios and mix thoroughly.

- 70% Ethanol - 7 parts 100% ethanol to 3 parts purified water
- 85% Ethanol - 8.5 parts 100% ethanol to 1.5 parts purified water

Store the solutions for up to 6 months at room temperature in an airtight container.

##### 2xSSC Solution

Dilute 1 part 20xSSC Solution with 9 parts purified water and mix thoroughly. Check pH and adjust to pH 7.0 using NaOH or HCl as required. Store the solution for up to 4 weeks at room temperature in an airtight container.

##### 0.4xSSC Solution

Dilute 1 part 20xSSC Solution with 49 parts purified water and mix thoroughly. Check pH and adjust to pH 7.0 using NaOH or HCl as required. Store the solution for up to 4 weeks at room temperature in an airtight container.

##### 2xSSC, 0.05% Tween-20 Solution

Dilute 1 part 20xSSC Solution with 9 parts purified water. Add 5µl of Tween-20 per 10ml and mix thoroughly. Check pH and adjust to pH 7.0 using NaOH or HCl as required. Store the solution for up to 4 weeks at room temperature in an airtight container.

#### FISH Protocol

(Note: Ensure that exposure of the probe and counterstain to laboratory lights is limited at all times).

##### Slide Preparation

- Spot the cell sample onto a glass microscope slide. Allow to dry. (**Optional, if using a cytogenetic drying chamber:** slides should be spotted using a cytogenetic drying chamber. The chamber should be operated at approximately 25°C and 50% humidity for optimal cell sample spotting. If a cytogenetic drying chamber is not available, use a fume hood as an alternative).
- Immerse the slide in 2xSSC for 2 minutes at room temperature (RT) without agitation.
- Dehydrate in an ethanol series (70%, 85% and 100%), each for 2 minutes at RT.
- Allow to dry.

##### Pre-Denaturation

- Remove the probe from the freezer and allow it to warm to RT. Briefly centrifuge tubes before use.
- Ensure that the probe solution is uniformly mixed with a pipette.
- Remove 10µl of probe per test, and transfer it to a microcentrifuge tube. Quickly return the remaining probe to the freezer.
- Place the probe and the sample slide to prewarm on a 37°C (+/- 1°C) hotplate for 5 minutes.
- Spot 10µl of probe mixture onto the cell sample and carefully apply a coverslip. Seal with rubber solution glue and allow the glue to dry completely.

##### Denaturation

- Denature the sample and probe simultaneously by heating the slide on a hotplate at 75°C (+/- 1°C) for 2 minutes.

##### Hybridisation

- Place the slide in a humid, lightproof container at 37°C (+/- 1°C) overnight.

##### Post-Hybridisation Washes

- Remove the DAPI from the freezer and allow it to warm to RT.
- Remove the coverslip and all traces of glue carefully.
- Immerse the slide in 0.4xSSC (pH 7.0) at 72°C (+/- 1°C) for 2 minutes without agitation.
- Drain the slide and immerse it in 2xSSC, 0.05% Tween-20 at RT (pH 7.0) for 30 seconds without agitation.
- Drain the slide and apply 10µl of DAPI antifade onto each sample.
- Cover with a coverslip, remove any bubbles and allow the colour to develop in the dark for 10 minutes.
- View with a fluorescence microscope. (See **Fluorescence Microscope Recommendation.**)

##### Stability of Finished Slides

Finished slides remain analysable for up to 1 month if stored in the dark at/or below RT.

##### Procedural Recommendations

- Baking or ageing of slides may reduce signal fluorescence
- Hybridisation conditions may be adversely affected by the use of reagents other than those provided or recommended by Cytocell Ltd
- Use a calibrated thermometer for measuring temperatures of solutions, waterbaths and incubators as these temperatures are critical for optimum product performance.
- The wash concentrations, pH and temperatures are important as low stringency can result in non-specific binding of the probe and too high stringency can result in a lack of signal
- Incomplete denaturation can result in lack of signal and over denaturation can also result in non-specific binding
- Over hybridisation can result in additional or unexpected signals
- Users should optimise the protocol for their own samples prior to using the test for diagnostic purposes

8. Suboptimal conditions may result in non-specific binding that may be misinterpreted as a probe signal

### Interpretation of Results





#### Assessing Slide Quality

The slide should not be analysed if:

- Signals are too weak to analyse in single filters - in order to proceed with analysis, signals should appear bright, distinct and easily evaluable
- There are high numbers of clumped/overlapping cells obstructing the analysis
- >50% of the cells are not hybridised
- There is excess of fluorescent particles between cells and/or a fluorescent haze that interferes with the signals - in optimal slides the background should appear dark or black and clean
- Cell nucleus borders cannot be distinguished and are not intact


#### Analysis Guidelines

- Two analysts should analyse and interpret each sample. Any discrepancies should be resolved by assessment by a third analyst
- Each analyst should be suitably qualified according to recognised national standards
- Each analyst should score independently 100 nuclei for each sample. The first analyst should start the analysis from the left side of the slide and the second analyst from the right one
- Each analyst should document their results in separate sheets
- Analyse only intact nuclei, not overlapped or crowded nuclei or nuclei covered by cytoplasmic debris or high degree of autofluorescence
- Avoid areas where there is excess of cytoplasmic debris or non-specific hybridisation
- Signal intensity may vary, even with a single nucleus. In such cases, use single filters and/or adjust the focal plane
- In suboptimal conditions signals may appear diffuse. If two signals of the same colour touch each other, or the distance between them is no greater than two signal widths, or when there is a faint strand connecting the two signals, count as one signal
- If in doubt about whether a cell is analysable or not, then do not analyse it

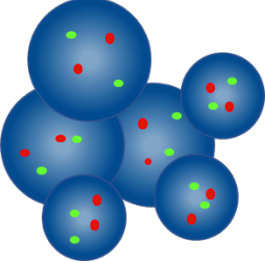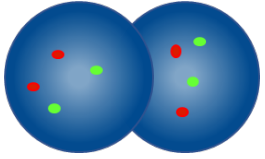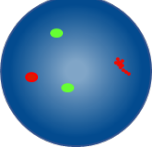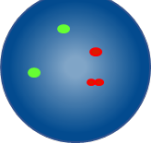Analysis Guidelines	
	Do not count – nuclei are too close together to determine boundaries
	Do not count overlapping nuclei – all areas of both nuclei are not visible
	Count as two red signals and two green signals – one of the two red signals is diffuse
	Count as two red signals and two green signals – the gap in one red signal is less than two signal widths


### Expected Results

#### Expected Normal Signal Pattern

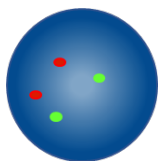

In a normal cell, two red and two green signals (2R, 2G) are expected.

#### Expected Abnormal Signal Pattern

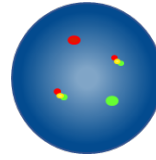

In a cell with a t(12;21)(p13.2;q22.1) translocation, the expected signal pattern will be one red, one green and two fusions (1R, 1G, 2F).

Other signal patterns are possible in aneuploid/unbalanced specimens.

#### Known Cross-Reactivity

No known cross-reactivity.

#### Adverse Event Reporting

If you believe this device has malfunctioned or suffered a deterioration in its performance characteristics which may have contributed to an adverse event (e.g. delayed or misdiagnosis, delayed or inappropriate treatment), this must be reported immediately to the manufacturer (**email:** [vigilance@ogt.com](mailto:vigilance@ogt.com)).

If applicable, the event should also be reported to your national competent authority. A list of vigilance contact points can be found at: <http://ec.europa.eu/growth/sectors/medical-devices/contacts/>.

#### Specific Performance Characteristics

##### Analytical Specificity

Analytical specificity is the percentage of signals that hybridise to the correct locus and no other location. The analytical specificity was established by analysing a total of 200 target loci. The analytical specificity was calculated as the number of FISH signals that hybridised to the correct locus divided by the total number of FISH signals hybridised.

Table 1. Analytical Specificity for the TEL/AML1 Translocation, Dual Fusion Probe

Probe	Target Locus	No. of Signals Hybridised to the Correct Locus	Total No. of Signals Hybridised	Specificity (%)
Red TEL	12p13.2	200	200	100
Green AML1	21q22.12	200	200	100

##### Analytical Sensitivity

Analytical sensitivity is the percentage of scoreable interphase cells with the expected normal signal pattern. The analytical sensitivity was established by analysing interphase cells across different normal samples. The sensitivity was calculated as the percentage of scoreable cells with the expected signal pattern (with a 95% confidence interval).

Table 2. Analytical Sensitivity for the TEL/AML1 Translocation, Dual Fusion Probe

No. of Cells with Expected Signal Patterns	No. of Cells with Scoreable Signals	Sensitivity (%)	95% Confidence Interval
475	499	95.0	1.6

##### Characterisation of Normal Cut-off Values

The normal cut-off value, in association with FISH probes, is the maximum percentage of scoreable interphase cells with a specific abnormal signal pattern at which a sample is considered normal for that signal pattern.

The normal cut-off value was established using samples from normal and positive patients. For each sample, the signal patterns of 100 cells were recorded. The Youden index was calculated to find the threshold value for which Sensitivity + Specificity-1 is maximised.

Table 3. Characterisation of Normal Cut-off Values for the TEL/AML1 Translocation, Dual Fusion Probe

Abnormal signal pattern	Youden Index	Normal Cut-off (%)
1R, 1G, 2F	1.00	3

Laboratories must verify cut-off values using their own data<sup>7, 8</sup>.

##### Precision and Reproducibility

Precision is a measure of the natural variation of a test when repeated several times under the same conditions. This was assessed by analysing repeats of the same lot number of probe tested on the same sample, in the same conditions on the same day.

Reproducibility is a measure of the variability of a test and has been established in terms of sample-to-sample, day-to-day and batch-to-batch variability. Day-to-day reproducibility was assessed by analysing the same samples on three different days. Batch-to-batch reproducibility was assessed by analysing the same samples using three different lot numbers of probe on one day. Sample-to-sample reproducibility was assessed by analysing three replicates of a sample on one day.

For each sample, signal patterns of 100 interphase cells were recorded and the percentage of cells with the expected signal pattern was calculated.

The reproducibility and precision were calculated as the Standard Deviation (STDEV) between replicates for each variable and overall mean STDEV.

**Table 4. Reproducibility and Precision for the TEL/AML1 Translocation, Dual Fusion Probe**

Variable	Standard Deviation (STDEV)
Precision	0.00
Sample-to-sample	0.00
Day-to-day	0.00
Batch-to-batch	0.00
Overall deviation	0.00

#### Clinical Performance

The clinical performance was established on a representative sample of the intended population for the product. For each sample, the signal patterns of  $\geq 100$  interphase cells were recorded. A normal/abnormal determination was made by comparing the percentage of cells with the specific abnormal signal pattern to the normal cut-off value. The results were then compared to the known status of the sample.

The results of the clinical data were analysed in order to produce sensitivity, specificity and cut off values using a one-dimensional approach.

**Table 5. Clinical Performance for the TEL/AML1 Translocation, Dual Fusion Probe**

Variable	Result
Clinical Sensitivity (true positive rate, TPR)	100%
Clinical Specificity (true negative rate, TNR)	100%
False Positive rate (FPR) = 1 – Specificity	0%

#### Additional Information

For additional product information please contact the Cytocell Technical Support Department.

T: +44 (0)1223 294048




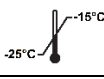


E: techsupport@cytozell.com

W: www.cytozell.com

#### References

- Mullighan, The Journal of Clinical Investigation 2012;122(12):3407-3415
- Wang *et al.*, Genes Dev 1998;12(15):2392-2402
- Swerdlow *et al.*, (eds.) WHO Classification of Tumours of Haematopoietic and Lymphoid Tissue, Lyon, France, 4th edition, IARC, 2017
- Borkhardt *et al.*, Blood. 1997;90(2):571-577
- Mosad *et al.*, Journal of Haematology & Oncology 2008;1:17
- Raynaud *et al.*, Blood 1996;87(7):2891-2899
- Arsham, MS., Barch, MJ. and Lawce HJ. (eds.) (2017) *The AGT Cytogenetics Laboratory Manual*. New Jersey: John Wiley & Sons Inc.
- Mascarello JT, Hirsch B, Kearney HM, et al. Section E9 of the American College of Medical Genetics technical standards and guidelines: fluorescence in situ hybridization. Genet Med. 2011;13(7):667-675.
- Wiktor AE, Dyke DLV, Stupca PJ, Ketterling RP, Thorland EC, Shearer BM, Fink SR, Stockero KJ, Majorowicz JR, Dewald GW. *Preclinical validation of fluorescence in situ hybridization assays for clinical practice*. Genetics in Medicine. 2006;8(1):16–23.

#### Guide to Symbols

REF	en: Catalogue number
IVD	en: <i>In vitro</i> diagnostic medical device
LOT	en: Batch code
	en: Consult instructions for use
	en: Manufacturer
	en: Use-by date
	en: Temperature limit
	en: Keep away from sunlight
	en: Contains sufficient for <n> tests
CONT	en: Contents

#### Patents and Trademarks

Aquarius and Cytocell are registered trademarks of Cytocell Ltd.



#### Cytocell Ltd.

Oxford Gene Technology,  
418 Cambridge Science Park,  
Milton Road,  
Cambridge, CB4 0PZ, UK  
T: +44(0)1223 294048  
F: +44(0)1223 294986  
E: probes@cytozell.com  
W: www.cytozell.com