



Instructions For Use

REF: LPH 079-S/LPH 079

E2A (TCF3)/PBX1 Translocation, Dual Fusion Probe



PROFESSIONAL USEONLY



Further information and other languages available at www.ogt.com

Limitations

This device is designed to detect rearrangements with breakpoints in the region bounded by the red and green clones in this probe set, which includes the *TCF3* and *PBX1* regions. Breakpoints outside this region or variant rearrangements wholly contained within this region may not be detected with this product.

The test is not intended for: use as a stand-alone diagnostic, prenatal testing, population-based screening, near-patient testing or self-testing. This product is intended for laboratory professional use only; all results should be interpreted by suitably-qualified staff, taking into account other relevant test results.

This product has not been validated for use on sample types or disease types other than those specified in the intended use.

Reporting and interpretation of FISH results should be consistent with professional standards of practice and should take into consideration other clinical and diagnostic information. This kit is intended as an adjunct to other diagnostic laboratory tests and therapeutic action should not be initiated on the basis of the FISH result alone.

Failure to adhere to the protocol may affect the performance and lead to false positive/negative results.

This kit has not been validated for purposes outside of the intended use stated.

Intended Use

The CytoCell E2A (TCF3)/PBX1 Translocation, Dual Fusion Probe is a qualitative, non-automated, fluorescence in situ hybridisation (FISH) test used to detect chromosomal rearrangements in the 1q23.3 region on chromosome 1 and 19p13.3 region on chromosome 19 in Carnoy's solution (3:1 methanol/acetic acid) fixed haematologically-derived cell suspensions from patients with confirmed or suspected acute lymphoblastic leukaemia (ALL).

Indications

This product is designed as an adjunct to other clinical and histopathological tests in recognised diagnostic and clinical care pathways, where knowledge of TCF3-PBX1 translocation status would be important for clinical management.

Principles of the test

Fluorescence in situ hybridisation (FISH) is a technique that allows DNA sequences to be detected on metaphase chromosomes or in interphase nuclei from fixed cytogenetic samples. The technique uses DNA probes that hybridise to entire chromosomes or single unique sequences, and serves as a powerful adjunct to Gbanded cytogenetic analysis. This technique can now be applied as an essential investigative tool within prenatal, haematological and solid tumour chromosomal analysis. Target DNA, after fixation and denaturation, is available for annealing to a similarly denatured, fluorescently labelled DNA probe, which has a complementary sequence. Following hybridisation, unbound and non-specifically bound DNA probe is removed and the DNA is counterstained for visualisation. Fluorescence microscopy then allows the visualisation of the hybridised probe on the target material.

Probe Information


The TCF3 (transcription factor 3) gene is located at 19p13.3 and PBX1 (PBX homeobox 1) at 1q23.3. Translocations involving TCF3 are some of the most common rearrangements in childhood B-cell acute lymphoblastic leukaemia (ALL)^{1,2}.

Two of the main TCF3 partners are PBX1 at 1q23.3 and HLF at 17q22. These become fused to TCF3 as a result of the t(1;19)(q23;p13) and t(17;19)(q22;p13) translocations, forming the TCF3-PBX1 and TCF3-HLF fusion genes, respectively. A rare cryptic inversion, inv(19)(p13;q13), has been reported to fuse TCF3 to TFPT (TCF3 fusion partner), resulting in the TCF3-TFPT fusion gene¹.

UK and European best practice guidelines suggest that when a TCF3 rearrangement is identified in B-cell ALL, it is important to distinguish between t(17;19)(q22;p13) and t(1;19)(q23;p13) as the former translocation is associated with adverse prognosis^{3,4}.

Probe Specification

PBX1, 1q23.3, Red E2A, 19p13.3, Green

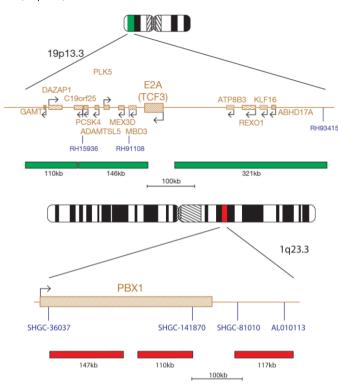

The E2A probe, labelled in green, contains two probes (110kb and 146kb) that cover the 3' end of the E2A (TCF3) gene and flanking region and a 321kb probe that covers a region 5' (centromeric) to the gene. The PBX1 probe, labelled in red, contains two probes (147kb and 110kb) that map within the PBX1 gene and a 117kb probe that maps 3' (telomeric) to the gene.

Materials Provided

Probe: 50µl per vial (5 tests) or 100µl per vial (10 tests)

The probes are provided premixed in hybridisation solution (formamide; dextran sulphate; saline-sodium citrate (SSC)) and are ready to use.

Counterstain: 150µl per vial (15 tests)

The counterstain is DAPI antifade (ES: 0.125µg/ml DAPI (4,6-diamidino-2-phenylindole)).

Warnings and Precautions

- For in vitro diagnostic use. For professional use only.
- 2. Wear gloves when handling DNA probes and DAPI counterstain.
- Probe mixtures contain formamide, which is a teratogen; do not breathe fumes
 or allow skin contact. Wear gloves, a lab coat, and handle in a fume hood.
 Upon disposal, flush with a large volume of water.
- DAPI is a potential carcinogen. Handle with care; wear gloves and a lab coat. Upon disposal, flush with a large volume of water.
- Dispose of all hazardous materials according to your institution's guidelines for hazardous waste disposal.
- 6. Operators must be capable of distinguishing the colours red, blue and green.
- Failure to adhere to the outlined protocol and reagents may affect the performance and lead to false positive/negative results.
- 8. The probe should not be diluted or mixed with other probes.
- Failure to use 10µl of probe during the pre-denaturation stage of the protocol
 may affect the performance and lead to false positive/negative results.

Storage and Handling



The kit should be stored between -25°C to -15°C in a freezer until the expiry date indicated on the kit label. The probe and counterstain vials must be stored in the dark.



The probe remains stable throughout the freeze-thaw cycles experienced during normal use (where one cycle constitutes the probe's removal from and replacement into the freezer) and is photostable for up to 48 hours after being exposed to continuous lighting conditions. All efforts must be made to limit exposure to light and temperature changes.

Equipment and Materials Necessary but not Supplied

Calibrated equipment must be used:

- Hotplate (with a solid plate and accurate temperature control up to 80°C)
- Calibrated variable volume micropipettes and tips range $1\mu I$ $200\mu I$
- Water bath with accurate temperature control at 37°C and 72°C
- Microcentrifuge tubes (0.5ml)
- Fluorescence microscope (Please see Fluorescence Microscope Recommendation section)
- Phase contrast microscope
- Clean plastic, ceramic or heat-resistant glass Coplin jars
- Forceps
- Calibrated pH meter (or pH indicator strips capable of measuring pH 6.5 -
- 10. Humidified container
- Fluores cence grade micros cope lens immersion oil
- Bench top centrifuge
- 13 Microscope slides
- 24x24mm coverslips
- 15. Timer
- 37°C incubator
- Rubber solution glue
- Vortex mixer 18.
- Graduated cylinders
- Magnetic stirrer
- 21. Calibrated thermometer

Optional Equipment not Supplied

Cytogenetic drying chambe

Reagents Needed but not Supplied

- 20x saline-sodium citrate (SSC) Solution
- 100% Ethanol
- Tween-20
- 1M Sodium hydroxide (NaOH)
- 1M Hydrochloric acid (HCI)
- Purified water

Fluorescence Microscope Recommendation

Use a 100-watt mercury lamp or equivalent and oil immersion plan apochromat objectives 60/63x or 100x for optimal visualisation. The fluorophores used in this probe set will excite and emit at the following wavelengths:

Fluorophore	Excitation _{max} [nm]	Emission _{max} [nm]
Green	495	521
Red	596	615

Ensure appropriate excitation and emission filters that cover the wavelengths listed above are fitted to the microscope. Use a triple bandpass DAPI/green spectrum/red spectrum filter or a dual bandpass green spectrum/red spectrum filter for optimal simultaneous visualisation of the green and red fluorophores.

Check the fluorescence microscope before use to ensure it is operating correctly. Use immersion oil that is suitable for fluorescence microscopy and formulated for low autofluorescence. Avoid mixing DAPI antifade with microscope immersion oil as this will obscure signals. Follow manufacturers' recommendations in regards to the life of the lamp and the age of the filters.

Sample Preparation

The kit is designed for use on haematologically-derived cell suspensions fixed in Carnoy's solution (3:1 methanol/acetic acid) fixative, that are prepared according to the laboratory or institution guidelines. Prepare air dried samples on microscope slides according to standard cytogenetic procedures. The AGT Cytogenetics Laboratory Manual contains recommendations for specimen collection, culturing, harvesting and for slide making⁵.

Solution Preparation

Ethanol Solutions

Dilute 100% ethanol with purified water using the following ratios and mix thoroughly:

- 70% Ethanol 7 parts 100% ethanol to 3 parts purified water
- 85% Ethanol 8.5 parts 100% ethand to 1.5 parts purified water

Store the solutions for up to 6 months at room temperature in an airtight container.

2xSSC Solution

Dilute 1 part 20xSSC Solution with 9 parts purified water and mix thoroughly. Check pH and adjust to pH 7.0 using NaOH or HCl as required. Store the solution for up to 4 weeks at room temperature in an airtight container.

0.4xSSC Solution

Dilute 1 part 20xSSC Solution with 49 parts purified water and mix thoroughly. Check pH and adjust to pH 7.0 using NaOH or HCl as required. Store the solution for up to 4 weeks at room temperature in an airtight container.

2xSSC, 0.05% Tween-20 Solution

Dilute 1 part 20xSSC Solution with 9 parts purified water. Add 5µl of Tween-20 per 10ml and mix thoroughly. Check pH and adjust to pH 7.0 using NaOH or HCl as required. Store the solution for up to 4 weeks at room temperature in an airtight container.

FISH Protocol

(Note: Ensure that exposure of the probe and counterstain to laboratory lights is limited at all times).

Slide Preparation

- Spot the cell sample onto a glass microscope slide. Allow to dry. (Optional, if using a cytogenetic drying chamber: slides should be spotted using a cytogenetic drying chamber. The chamber should be operated at approximately 25°C and 50% humidity for optimal cell sample spotting. If a cytogenetic drying chamber is not available, use a fume hood as an alternative).
- Immerse the slide in 2xSSC for 2 minutes at room temperature (RT) without agitation.
- 3. Dehydrate in an ethanol series (70%, 85% and 100%), each for 2 minutes at RT.
- Allow to dry.

Pre-Denaturation

- Remove the probe from the freezer and allow it to warm to room temperature (RT). Briefly centrifuge tubes before use.
- Ensure that the probe solution is uniformly mixed with a pipette.
- Remove 10µl of probe per test, and transfer it to a microcentrifuge tube. Quickly return the remaining probe to the freezer.
- Place the probe and the sample slide to prewarm on a 37°C (+/- 1°C) hotplate for 5 minutes.
- Spot 10µl of probe mixture onto the cell sample and carefully apply a coverslip. Seal with rubber solution glue and allow the glue to dry completely.

10. Denature the sample and probe simultaneously by heating the slide on a hotplate at 75°C (+/- 1°C) for 2 minutes.

Hybridisation

11. Place the slide in a humid, lightproof container at 37°C (+/- 1°C) overnight.

Post-Hybridisation Washes

- 12. Remove the DAPI from the freezer and allow it to warm to room temperature
- Remove the coverslip and all traces of glue carefully.
- Immerse the slide in 0.4xSSC (pH 7.0) at 72°C (+/- 1°C) for 2 minutes without agitation.
- Drain the slide and immerse it in 2xSSC, 0.05% Tween-20 at RT (pH 7.0) for 30 seconds without agitation.
- Drain the slide and apply 10µl of DAPI antifade onto each sample.
- Cover with a coverslip, remove any bubbles and allow the colour to develop in the dark for 10 minutes.
- View with a fluorescence microscope. (See Fluorescence Microscope Recommendation.)

Stability of Finished Slides

Finished slides remain analysable for up to 1 month if stored in the dark at/or below

Procedural Recommendations

- Baking or ageing of slides may reduce signal fluorescence Hybridisation conditions may be adversely affected by the use of reagents other than those provided or recommended by Cytocell Ltd
- Use a calibrated thermometer for measuring temperatures of solutions, waterbaths and incubators as these temperatures are critical for optimum product performance.
- The wash concentrations, pH and temperatures are important as low stringency can result in non-specific binding of the probe and too high stringency can result in a lack of signal
- Incomplete denaturation can result in lack of signal and over denaturation can also result in non-specific binding
- Over hybridisation can result in additional or unexpected signals
- Users should optimise the protocol for their own samples prior to using the test for diagnostic purposes
- Suboptimal conditions may result in non-specific binding that may be misinterpretedas a probe signal


Interpretation of Results

Assessing Slide Quality

- The slide should not be analysed if: Signals are too weak to analyse in single filters - in order to proceed with analysis, signals should appear bright, distinct and easily evaluable
- There are high numbers of clumped/overlapping cells obstructing the analysis
- >50% of the cells are not hybridised
- There is excess of fluorescent particles between cells and/or a fluorescent haze that interferes with the signals - in optimal slides the background should appear dark or black and clean
- Cell nucleus borders cannot be distinguished and are not intact

Analysis Guidelines

- Two analysts should analyse and interpret each sample. Any discrepancies should be resolved by assessment by a third analyst
- Each analyst should be suitably qualified according to recognised national standards
- Each analyst should score independently 100 nuclei for each sample. The first analyst should start the analysis from the left side of the slide and the second analyst from the right one
- Each analyst should document their results in separate sheets
- Analyse only intact nuclei, not overlapped or crowded nuclei or nuclei covered by cytoplasmic debris or high degree of autofluorescence
- Avoid areas where there is excess of cytoplasmic debris or non-specific hybridisation
- Signal intensity may vary, even with a single nucleus. In such cases, use single filters and/or adjust the focal plane
- In suboptimal conditions signals may appear diffuse. If two signals of the same colour touch each other, or the distance between them is no greater than two signal widths, or when there is a faint strand connecting the two signals, count as one signal
- If in doubt about whether a cell is analysable or not, then do not analyse it

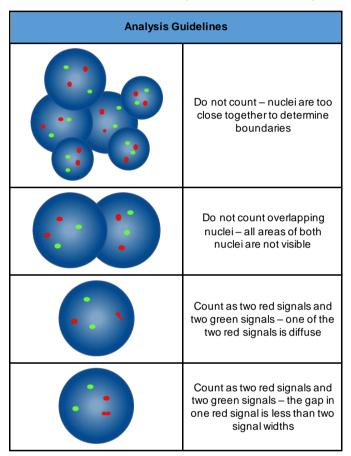

Expected Results
Expected Normal Signal Pattern



In a normal cell, two red and two green signals (2R, 2G) are expected.

Expected Abnormal Signal Patterns



In a cell with a balanced t(1;19)(q23;p13.3), the expected signal pattern will be one red, one green and two fusion signals which may appear as red and green signals close together (1R, 1G, 2F).



In a cell with an unbalanced der(19)t(1;19)(q23;p13.3), the expected signal pattern will be two red, one green and one fusion signal (2R, 1G, 1F).

Other signal patterns are possible in aneuploid/unbalanced specimens.

Known Cross-Reactivity

No known cross-reactivity.

Adverse Event Reporting

If you believe this device has malfunctioned or suffered a deterioration in its performance characteristics which may have contributed to an adverse event (e.g. delayed or misdiagnosis, delayed or inappropriate treatment), this must be reported immediately to the manufacturer (email: vigilance@ogt.com).

If applicable, the event should also be reported to your national competent authority. A list of vigilance contact points can be found at: http://ec.europa.eu/growth/sectors/medical-devices/contacts/.

Specific Performance Characteristics Analytical Specificity

Analytical specificity is the percentage of signals that hybridise to the correct locus and no other location. The analytical specificity was established by analysing a total of 200 target loci. The analytical specificity was calculated as the number of FISH signals that hybridised to the correct locus divided by the total number of FISH signals hybridised.

Table 1. Analytical Specificity for the E2A/PBX1 Translocation, Dual Fusion Probe

Probe	Target Locus	No. of Signals Hybridised to the Correct Locus	Total No. of Signals Hybridised	Specificity (%)
E2A Green	19p13.3	200	200	100
PBX1 Red	1q23.3	200	200	100

Analytical Sensitivity

Analytical sensitivity is the percentage of scoreable interphase cells with the expected normal signal pattern. The analytical sensitivity was established by analysing interphase cells across different normal samples. The sensitivity was calculated as the percentage of scoreable cells with the expected signal pattern (with a 95% confidence interval).

Table 2. Analytical Sensitivity for the E2A/PBX1 Translocation, Dual Fusion Probe

No. of Cells with Expected Signal Patterns	No. of Cells with Scoreable Signals	Sensitivity (%)	95% Confidence Interval
478	500	95.6	1.5

Characterisation of Normal Cut-off Values

The normal cut-off value, in association with FISH probes, is the maximum percentage of scoreable interphase cells with a specific abnormal signal pattern at which a sample is considered normal for that signal pattern.

The normal cut-off value was established using samples from normal and positive patients. For each sample, the signal patterns of 100 cells were recorded. The Youden index was calculated to find the threshold value for which Sensitivity + Specificity-1 is maximised.

<u>Table 3. Characterisation of Normal Cut-off Values for the E2A/PBX1 Translocation, Dual Fusion Probe</u>

Abnormal signal pattern	Youden Index	Normal Cut-off (%)
2R, 1G, 1F	0.99	1

Laboratories must verify cut-off values using their own data^{6,7}.

Precision and Reproducibility

Precision is a measure of the natural variation of a test when repeated several times under the same conditions. This was assessed by analysing repeats of the same lot number of probe tested on the same sample, in the same conditions on the same day.

Reproducibility is a measure of the variability of a test and has been established in terms of sample-to-sample, day-to-day and batch-to-batch variability. Day-to-day reproducibility was assessed by analysing the same samples on three different days. Batch-to-batch reproducibility was assessed by analysing the same samples using three different lot numbers of probe on one day. Sample-to-sample reproducibility was assessed by analysing three replicates of a sample on one day. For each sample, signal patterns of 100 interphase cells were recorded and the percentage of cells with the expected signal pattern was calculated.

The reproducibility and precision were calculated as the Standard Deviation (STDEV) between replicates for each variable and overall mean STDEV.

Table 4. Reproducibility and Precision for the E2A/PBX1 Translocation, Dual Fusion Probe

Variable	Standard Deviation (STDEV)
Precision	0.00
Sample-to-sample	0.00
Day-to-day	0.00
Batch-to-batch	0.00
Overall deviation	0.00

Clinical Performance

The clinical performance was established on a representative sample of the intended population for the product. For each sample, the signal patterns of ≥100 interphase cells were recorded. A normal/abnormal determination was made by comparing the percentage of cells with the specific abnormal signal pattern to the normal cut-off value. The results were then compared to the known status of the sample.

The results of the clinical data were analysed in order to produce sensitivity, specificity and cut off values using a one-dimensional approach.

Table 5. Clinical Performance for the E2A/PBX1 Translocation, Dual Fusion Probe

Variable	Result
Clinical Sensitivity (true positive rate, TPR)	99.9%
Clinical Specificity (true negative rate, TNR)	100%
False Positive rate (FPR) = 1 – Specificity	0.00%

Additional Information

For additional product information please contact the CytoCell Technical Support Department.

T: +44 (0)1223 294048

E: techsupport@cytocell.com

W: www.ogt.com

References

- Van der Burg et al., Leukemia 2004;18(5):895-908
- Swerdlow et al., (eds,) WHO Classification of Tumours of Haematopoietic and Lymphoid Tissue, Lyon, France, 4th edition, IARC, 2017
- 3. Professional Guidelines for Clinical Cytogenetics: Acute Lymphoblastic Leukaemia BEST PRACTICE GUIDELINES (2011) www.cytogenetics.org.uk
- Hasting et al., Guidelines and Quality Assurance for Acquired Cytogenetics 4. (2013) ECA Newsletter 31
- Arsham, MS., Barch, MJ. and Lawce HJ. (eds.) (2017) The AGT Cytogenetics Laboratory Manual. New Jersey: John Wiley & Sons Inc. Mascarello JT, Hirsch B, Kearney HM, et al. Section E9 of the American
- College of Medical Genetics technical standards and guidelines: fluorescence in situ hybridization. Genet Med. 2011;13(7):667-675.
- Wiktor AE, Dyke DLV, Stupca PJ, Ketterling RP, Thorland EC, Shearer BM, Fink SR. Stockero KJ. Maiorowicz JR. Dewald GW. Preclinical validation of fluorescence in situ hybridization assays for clinical practice. Genetics in Medicine. 2006;8(1):16-23.

Guide to Symbols

REF	en: Catalogue number
IVD	en: In vitro diagnostic medical device
LOT	en: Batch code
[]i	en: Consult instructions for use
	en: Manufacturer
\square	en: Use-by date
-25°C	en: Temperature limit
茶	en: Keep away from sunlight
Σ	en: Contains sufficient for <n> tests</n>
CONT	en: Contents

Patents and Trademarks

CytoCell is a registered trademark of Cytocell Ltd.



Cytocell Ltd.
Oxford Gene Technology,
418 Cambridge Science Park, Milton Road, Cambridge, CB4 0PZ, UK T: +44(0)1223 294048 F: +44(0)1223 294986 E: probes @cytocell.com W: www.ogt.com