



A Sysmex Group Company



AML1 (RUNX1) Distal Probe Green


REF: LPH 573-A

Analyte Specific Reagent: Analytical and performance characteristics are not established.

Fluorescence *In Situ* Hybridisation (FISH) is a technique that allows DNA sequences to be detected on metaphase chromosomes or in interphase nuclei from fixed cytogenetic samples. The technique uses DNA probes that hybridise to entire chromosomes or single unique sequences, and serves as a powerful adjunct to classic cytogenetics. Recent developments have meant that this valuable technique can now be applied as an essential tool in prenatal, haematological and pathological chromosomal analysis. Target DNA, after fixation and denaturation, is available for annealing to a similarly denatured, fluorescently labelled DNA probe, which has a complementary sequence. Following hybridisation, unbound and non-specifically bound DNA probe is removed and the DNA is counterstained for visualisation. Fluorescence microscopy then allows the visualisation of the hybridised probe on the target material.

Probe Specification

AML1 (RUNX1) Distal, 21q22.1

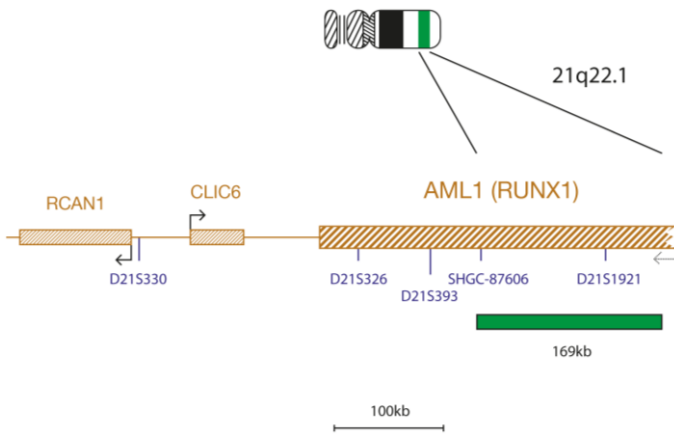

The CytoCell AML1 (RUNX1) Distal probe consists of a 169kb probe, labelled in green, covering part of the *AML1 (RUNX1)* gene, including markers SHGC-87606 and D21S1921.

Materials Provided

Probe: 50µl per vial

Probe concentration: 11-16.4ng/µl

The probe is provided in hybridisation solution (Formamide; Dextran Sulphate; SSC) and is ready to use.

Warnings and Precautions

1. For professional use only.
2. Wear gloves when handling DNA probes.
3. Probe contains formamide, which is a teratogen; do not breathe fumes or allow skin contact. Wear gloves, a lab coat, and handle in a fume hood. Upon disposal, flush with a large volume of water.
4. Dispose of all hazardous materials according to your institution's guidelines for hazardous waste disposal.
5. Operators must be capable of visually distinguishing between red, blue and green.

Storage and Handling

Store the probe at -20°C until the expiry date indicated on the label. Store the probe vial in the dark. Ensure that exposure of the probe to laboratory lights is limited at all times.

Known Cross-Reactivity

No known cross-reactivity.

Additional Information

For additional product information please contact the CytoCell Technical Support Department.

T: +44 (0)1223 294048

E: techsupport@cytoCell.com

W: www.ogt.com

Patents and Trademarks

CytoCell is a registered trademark of CytoCell Ltd.



CytoCell Ltd.

Oxford Gene Technology,
418 Cambridge Science Park,
Milton Road,

Cambridge, CB4 0PZ, UK

T: +44(0)1223 294048

F: +44(0)1223 294986

E: probes@cytoCell.com

W: www.ogt.com